

ASTM D2866

Milestone Application Note for Ashing – UltraFAST - Officials



Summary

This procedure covers the determination of total ash content of activated carbon samples according to ASTM D2866 method.

Instrumentation

Milestone PYRO microwave system


UltraFAST muffle furnace

Crucibles (ceramic, quartz, quartz fibre)

Analytical balance ($\pm 0,0001\text{g}$)

Disseccator

Spatula

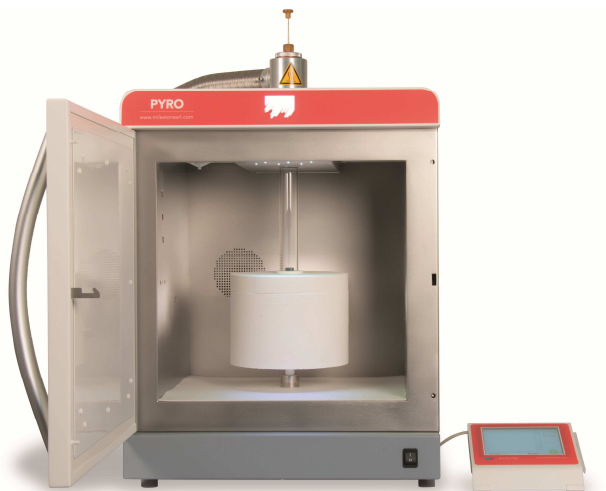

Sample amount

Sample weight	1-5 g (according to the expected ash content)
---------------	---

Procedure

1. Place a clean crucible on the balance plate, tare it and weigh the sample.
2. Place the crucible into the preheated furnace.
3. Remove the crucible from the furnace and cool it in dissector.
4. Place the crucible on the balance plate and take the weigh.
5. Calculate the ash residue.

Notes

The crucible before use must be at constant weight.

All modifications to the present method must be agreed with Milestone Application Specialist.

General precautions

Always use hand, eye and body protection when operating with the microwave system.

Microwave program

Step	Time	T1	Power	Vent.
1	00:06:00	650°C	1800 W	*
2	02:00:00*	650°C	1800 W	*

* Ashing time must be extended until all the carbon is removed.

Subject to change without notice.
For additional information please contact
application@milestonesrl.com

